



PRO-DENSE™ BONE GRAFT SUBSTITUTE

MIXING INSTRUCTIONS

154101-1


The following languages are included in this packet:

English (en)	Deutsch (de)	Nederlands (nl)	Français (fr)
Español (es)	Italiano (it)	Português (pt)	Türkçe (tk)

For additional languages, visit our website www.wright.com.

Then click on the **Prescribing Information** option.

For additional information and translations please contact the manufacturer or local distributor.

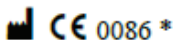

Wright Medical Technology, Inc.
1023 Cherry Road
Memphis, TN 38117
U.S.A.

EC	REP
----	-----

Tornier S.A.S
161 Rue Lavoisier
38330 Montbonnot Saint Martin
France

*** The CE-Marking of Conformity is applied per catalog number and appears on the outer label, if applicable.**

R_x ONLY

July 2018
Printed in U.S.A.

NOTE: Do not add any additional substances to the graft. Using alternative mixing solutions and/or adding other substances to the mixture may alter the safety, effectiveness and setting of this product.

1. Remove metal crimp from solution vial. Attach mixing syringe to powder container with **RED** adaptor and pull back the mixing rod.




2. Invert container and tap to transfer ALL powder into mixing syringe. Remove powder container.




3. Remove the vial stopper and add ALL solution into mixing syringe.



4. Remove RED adaptor and attach WHITE mixing cap. Remove the plunger sleeve by disconnecting the handle end from the mixing rod first.

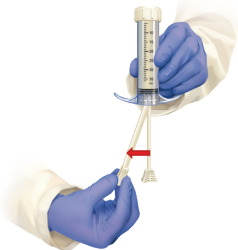

5. Mix for 30 seconds by moving the mixing rod back and forth while rotating. **Be sure to break-up and incorporate the powder 'plug' near the back and front of the syringe.**

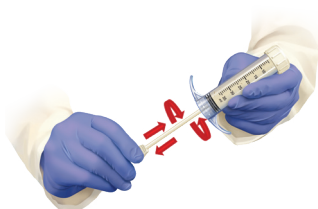

Note: Inject paste within 5 min.

6. **Completely** pull back the mixing rod and remove the mixing cap. Reattach the plunger sleeve onto the mixing rod by positioning the front end onto the piston first before snapping into place.



7. Connect delivery syringe to mixing syringe with **ORANGE** adaptor. Delivery syringe plunger should be **fully depressed**.



8. Gently push the mixing syringe plunger to transfer the graft into the delivery syringe. Trapped air will be removed during transfer. Repeat steps 7-8 for additional syringes.




9. Disconnect delivery syringe from the **ORANGE** adaptor and attach the provided luer cap and accessory needle. Inject graft into bone defect.



Note: The graft will set in ~ 20 - 30 min.

10. (OPTIONAL) For mechanical assist, remove plunger sleeve and snap off tip of mixing rod.

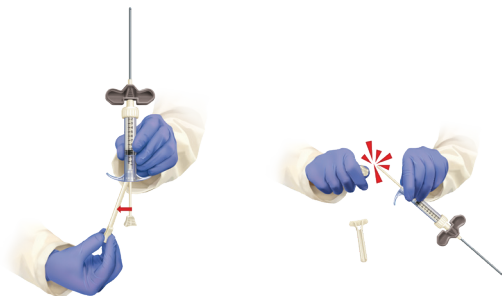

11. (OPTIONAL) Insert spindle drive over the mixing rod until you hear the spindle nut engage with both tabs on the syringe. Rotate handle clockwise to manually eject graft.



NOTE: Residual material may be left in syringe. The amount of remaining material is relative to user's adherence to the mixing and ejecting times outlined in these instructions.